

Recombinant tuna FGF-2 (154 aa) protein (Qk105)



Type: Stem cells

Available for purchase: Unit Size (µg): 25, 50, 100, 500, 1000

Buy online with secure credit card or purchase order.

[View this product and buy online](#)

Product Information

Recombinant tuna FGF-2 protein 154 aa (bFGF/basic FGF) for the development of optimized serum-free culture media for species-specific yellowfin tuna (*Thunnus albacares*) in [cellular agriculture](#) protocols and veterinary research applications. This longer form of FGF-2 is used in comparative cell culture media optimization studies alongside [Qk104](#), the 145 aa form of tuna FGF-2. FGF-2 is used extensively in the maintenance and proliferation of induced pluripotent (iPSC) and embryonic stem cells (ESC) and for enhancement of proliferation in primary tuna cell culture. Receptor binding affinity and efficacy may differ depending on each species. Using a species-specific growth factor enhances receptor binding affinity, resulting in a lower concentration required in culture.

Tuna FGF-2 is a high purity 17 kDa FGF-2 / bFGF protein, [animal origin-free](#) (AOF) and carrier-protein free (CF).

Alternative protein names

Basic fibroblast growth factor, bFGF, FGF-β, FGF2, FGF 2, Fibroblast growth factor-basic, HBGF-2, betaFGF, beta FGF

Molecular weight

17 kDa (monomer)

Protein Uniprot number

High purity tuna protein (Uniprot number: XP_044218972.1)

Species reactivity

- tuna

Product Information

- >98%, by SDS-PAGE quantitative densitometry
- Expressed in *E. coli*
- Animal origin-free (AOF) and carrier protein-free
- Manufactured in our Cambridge, UK laboratories
- Lyophilized from Tris, NaCl, Cys, mannitol

Reconstitution instructions

- Resuspend in sterile-filtered water at >50 µg/ml

Featured applications

- Expansion of tuna pluripotent, embryonic and mesenchymal stem cells
- Cellular agriculture and cultivated meat cell culture media optimization
- Serum-free media development
- Cellular proliferation, migration and survival


Further quality assays

- Mass spectrometry: single species with expected mass
- Recovery from stock vial: >95%
- Endotoxin: <0.05 EU/µg protein

Scientific Information


Bioactivity

Recombinant tuna FGF-2 154 aa activity was determined using the Promega serum response element luciferase reporter assay in transfected HEK293T cells. Cells were treated in triplicate with a serial dilution of FGF-2 for 3 hours. Firefly luciferase activity was measured and normalized to the control Renilla luciferase activity. Data from Qk105 lot 204640. EC50 = 0.139 ng/ml (8.2 pM).

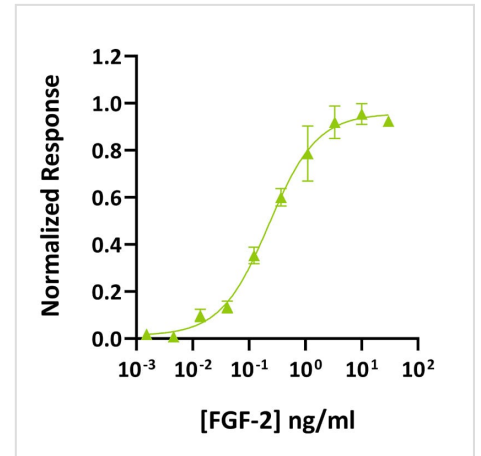

Purity

Recombinant tuna FGF-2 154 aa migrates as a major band at approximately 17 kDa (monomer) in reduced (R) and non-reduced (NR) conditions. No contaminating protein bands are present. The purified recombinant protein (3 µg) was resolved using 15% w/v SDS-PAGE in reduced (+β-mercaptoethanol, R) and non-reduced (NR) conditions and stained with Coomassie Brilliant Blue R250. Data from Qk105 lot #204640.

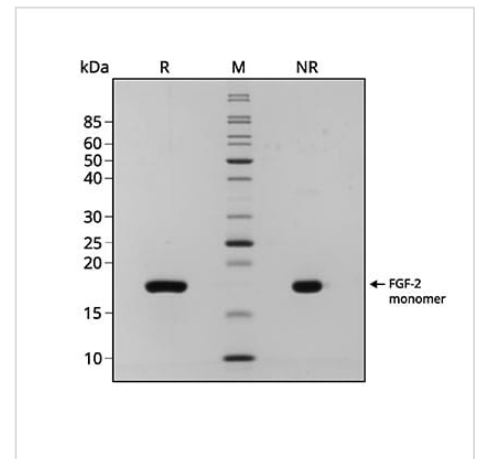

Original product page: <https://ryan.calliope-alpha.ts.net/product/recombinant-tuna-fgf-2-154aa-protein-qk105/>

PDF generated: 12 May 2026

Copyright © 2026 by Qkine Ltd. All rights reserved including graphics and images.